


PROGRESSIVE PureLife FREE® 2.

The best entry level progressive lens with unique benefits.


ZERO BASE CURVE EFFECT.

In comparison to other lenses in the branded entry segment or conventional lenses Rodenstock Progressive PureLife Free® eliminates the base curve effect through spherical optimization. Thus ever spherical prescription results in optimum optical performance for 100% satisfied customers.

Conventional / Standard

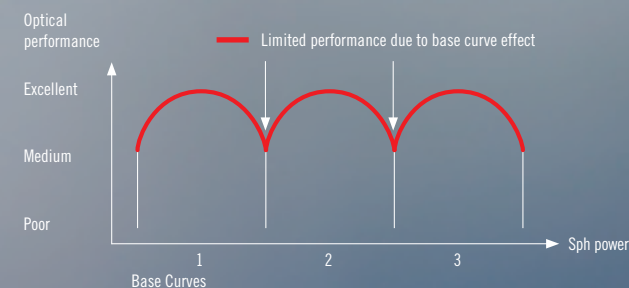

Progressiv PureLife Free® 2

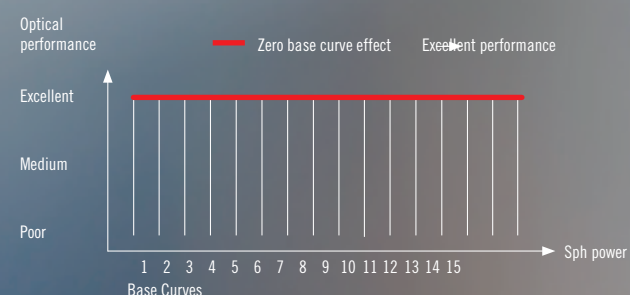

LENS DATA:

Progressiv	Progressiv Pure Life Free® 2 L	Progressiv Pure Life Free® 2 M	Progressiv Pure Life Free® 2 XS
Brand engraving	R	R	R
Functional engraving	○	○	○
Additional engraving	PL2	PM2	PX2
Progression zone length	18 mm	16 mm	14 mm
Glazing height	20 mm	18 mm	16 mm
Total height	28 mm	26 mm	24 mm
Base Curve/Index	under functional engraving nasal		
Addition	under functional engraving temporal		
Necessary order parameters			
Refraction data	far value and additon		
Optional order parameters/product data			
Frame-/Centering data	✓	✓	✓
Technical details			
<ul style="list-style-type: none"> For all lenses with Position of Wear: order values are printed on the envelope Adaption: straight forward position with a habitual head and body position 			

TECHNICAL DETAILS:

	Ranges from	Ranges to	cyl	Addition
Progressiv PureLife Free® 2 1.74 L/M/XS	-10.00	+8.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.67 ColorMatic IQ L/M/XS	-10.00	+8.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.67 L/M/XS	-10.00	+8.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.60 ColorMatic IQ L/M/XS	-10.00	+8.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.60 ColorMatic IQ Sun L/M/XS	-10.00	+8.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.6 Polarized L/M/XS	-6.00	+6.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.60 L/M/XS	-10.00	+8.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.59 Polycarbonate L/M/XS	-8.00	+4.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.54 ColorMatic IQ L/M/XS	-8.00	+6.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.53 Trivex L/M/XS	-8.00	+4.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.50 Polarized L/M/XS	-6.00	+6.00	+4.00	0.75-3.50
Progressiv PureLife Free® 2 1.50 L/M/XS	-8.00	+6.00	+4.00	0.75-3.50

All ranges show at a diameter of max. 70 mm. Deviation causes a change of +/- 2 D.

MATERIAL FACTS:

	Refractive index n _e	Abbè No., v _e	Density ρ (g/cm ²)	100% UV protection up to*	Tints
1.74	1.737	32.5	1.47	400	
1.67 CM IQ	1.668	31.4	1.37	400	Brown, grey 8-85%
1.67	1.668	31.4	1.37	400	
1.60 CM IQ	1.597	40.5	1.30	400	Brown, grey, green 8 - 85%
1.60 CM IQ Sun	1.597	40.5	1.30	400	All colours 55 - 90%
1.6 Polarized	1.6	42.0	1.30	400	Brown 78%, grey 83%, green 85%
1.60	1.597	40.5	1.30	400	
1.59 Polycarbonate	1.591	29.8	1.20	385	
1.54 CM IQ	1.540	43.0	1.19	400	Brown, grey 5-85%
1.53 Trivex	1.530	45.0	1.11	400	
1.50 Polarized	1.502	58.0	1.32	380	Brown 78%, grey 83%, green 85%
1.50	1.502	58.2	1.32	350	

*Data according to EN ISO 14889:2013 and 8980-3:2013 or 12312-1:2013 for an untreated lens without any coating.



PROGRESSIVE PureLife FREE® 2.

RODENSTOCK
SUPERIOR

The best entry level progressive lens with unique benefits

Rodenstock Progressiv PureLife Free® 2 is the true progressive lens for starters who want to rely on brand quality.

BENEFITS:

- With EyeModel up to 25% better vision in the near and intermediate range
- The most modern position of wear for better vision thanks to the Retina Focus Principle
- Eliminator of the base curve effect and very finely stepped base curve system for an aesthetically perfect fit
- Rodenstock quality standard
- Brilliant, clear vision from far to near
- Remarkable larger and more stabilized vision fields
- Fatigue-free reading
- Very good spontaneous compatibility
- Exceptional price/performance ratio and security of a known brand

FEATURES:

- EyeLT EyeModel
- Zero base curve effect
- Variable and power-dependent inset
- Finely stepped base curve system for an aesthetically perfect fit
- Retina-focus principle
- 3 D Freeform Technology

