



# A COMPLETE **SINGLE VISION** DESIGN FAMILY

 **DIGITAL** RAY-PATH® 2

 IOT  
**Endless**®


 IOT  
**Endless**® Plus

Single vision and anti-fatigue lenses utilizing **advanced technology**  
**for better vision** in dress, sun and sports eyewear.

  
**PLASTIC PLUS**  
*NO ONE DOES DIGITAL BETTER*

**Benefits**

- 1 Precise and comfortable focus at all working distances in any direction of gaze.
- 2 Virtual elimination of peripheral blur.
- 3 Customized for the individual parameters of each wearer, material, base curve, and frame.
- 4 Automatic centering for thickness reduction.

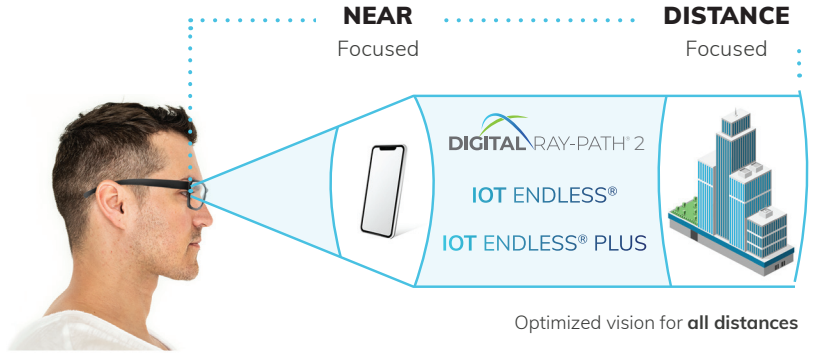

**Digital Ray-Path® 2 Technology personalizes your patient's lenses for crisp vision at all distances.** They will enjoy clearer vision with IOT Endless® and IOT Endless® Plus lenses.

**Impeccable visual quality at any distance**

IOT Endless® and IOT Endless® Plus lenses include Digital Ray-Path® 2 Technology. In addition to mathematically compensating for oblique aberrations, Digital Ray-Path® 2 adds the intelligent use of the wearer's accommodation; the small adjustments the eyes naturally make to view objects at different distances.

The result is a **drastic reduction of oblique aberrations** throughout the visual field. Wearers will enjoy **impeccable visual quality, greater comfort, and more precise focus.**



**Lenses for any single vision wearer**

<p><b>IOT Endless®</b></p> <p><b>0.00 D</b></p> <p><b>Ideal</b> wearer</p> <p>Single vision lens wearers of any age with <b>no symptoms of presbyopia or digital eyestrain.</b></p>	<p><b>IOT Endless® Plus</b></p> <p><b>0.50 D</b></p> <p><b>Ideal</b> wearer</p> <p>Single vision lens wearers with <b>mild symptoms of digital eyestrain.</b></p>	<p><b>IOT Endless® Plus</b></p> <p><b>0.75 D</b></p> <p><b>Ideal</b> wearer</p> <p>Single vision lens wearers with <b>moderate symptoms of digital eyestrain.</b></p>	<p><b>IOT Endless® Plus</b></p> <p><b>1.00 D</b></p> <p><b>Ideal</b> wearer</p> <p>Pre-presbyopes with <b>symptoms of digital eyestrain and difficulty focusing at near distances.</b></p>
---	---	---	--

IOT Endless® and IOT Endless® Plus lenses are fully personalized and compensated. If full personalization parameters are not provided, defaults of wrap 5°, panto 8°, and vertex 14 mm will be used. IOT Endless® and IOT Endless® Plus lenses are dispensed and fitted like traditional single vision lenses, however it is recommended to have a minimum fitting height of 16 mm for IOT Endless® Plus lenses. Contact Plastic Plus for more information.